


# DentuSil™ Silicone Soft Reline

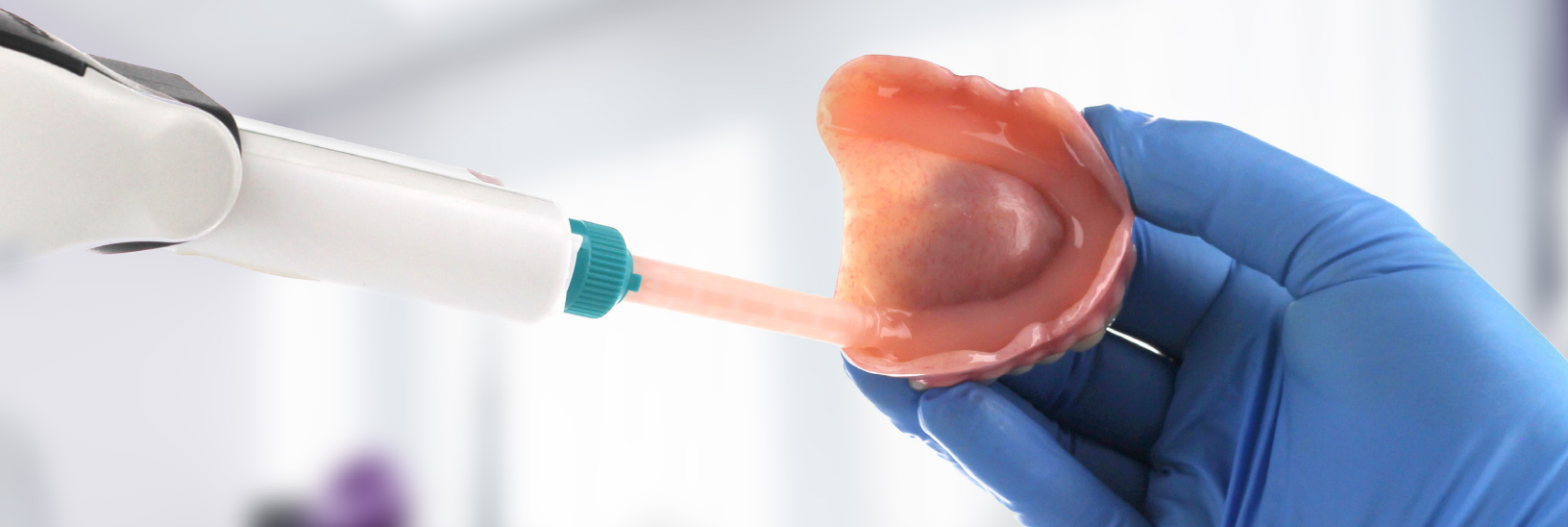

**DentuSil is a cartridge-dispensed silicone soft reline material that provides patients with the comfortable, reliable reline in one office visit.**

## About DentuSil

DentuSil's silicone base absorbs less water, making the reline more resistant to stains and odors. Its translucent pink color provides a natural, healthy appearance. The Adhesive is specially formulated to hold material in place, providing a secure reline in minutes.

## Technical Tips

Addition Reaction Silicone

Mixing Time: Automix

Working Time: 1 minute

Setting Time: 5-7 minutes





Polymerization Shrinkage: Less than 0.1%

Shelf Life: 2 years stored at or below 72°F (21°C)

## Features & Benefits

- Translucent pink color gives a natural, healthy appearance
- Low porosity material resisting stains and odors
- No monomer making it perfect for use on patients with sensitivity to acrylics

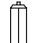

PRODUCT	COMPONENTS	ITEM NUMBER
Standard Kit	 50 ml Cartridge  Mixing Tips  10 ml Liquid Adhesive  10 ml Silicone Release	0921276
Silicone Release	10 ml	0921279
Applications	10 applications	
Shelf Life	2 years	
Mouth Life	1 year	
Shore A Hardness	25	



[dental.keystoneindustries.com](http://dental.keystoneindustries.com)

US: 480 South Democrat Road Gibbstown NJ, 08027

EU: Keystone Industries GmbH Stockholzstr. 11 78224 Singen, Germany

US: 1.800.333.3131  
EU: +49 77 31 91 21 01